
CE

⚠WARNING: Please mind potential fire risk and strictly follow the instruction manual for charging and storage. For disposal, please check your local authority's website for more information and dispose of the lithium-ion battery safely in accordance with the applicable legal requirements.

POWERED
EAR WAX REMOVER

USER MANUAL

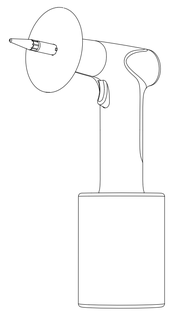Thank you for purchasing our product. In order to ensure safety and better perfmance,please carefully read this manual before use

Catalog

English 01-05
Deutsch 06-10
Français 11-15
Español 16-20
Italiano 21-25
日本語 26-30

CE FC UK CA RoHS Made in China

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Before use, please read this instruction manual carefully and completely, and keep it for future reference!

Precautions

Please read all instructions, warnings and tips carefully before use. Improper use may cause serious personal injury and property damage.


1.Do not put corrosive liquids into the product to avoid damage to the human body and the product.
2. It is normal to feel pressure but no pain. If you feel pain, please stop using this product immediately and consult your doctor for professional help.
3. Children or people with disabilities should use this product under the supervision of an adult or guardian.
4. This product is not a toy. Please keep it out of reach of children.
5. This product is equipped with a battery. Please keep it away from fire, moisture and high temperature. It may explode due to overheating.
6. If you do not use this product, please store it in a dry place.
7. Make sure to connect a suitable power supply DC5V to charge the product.
8. Do not use solutions below body temperature to avoid discomfort. It is recommended to use warm water (not too hot to avoid burns)
9. Do not change or modify the product, or disassemble the internal equipment, otherwise it may damage the product or cause personal injury.
10. Do not classify it as household waste. When you no longer need this product, please use your local government's recycling process and take it to a local electrical equipment recycling station.

Warm Tips: All pictures in the manual are for reference only. Due to product updates or upgrades, there may be differences between the actual product and the picture, please refer to the actual product.

Troubleshooting

1. If the product cannot be turned on, please charge it to prevent low battery.
2. If the product has a short working time after being fully charged, please charge and use it in an ambient temperature between 15~35℃, or if the charging time is insufficient, please charge it.
3. If the LED light does not light up when charging, please check whether the power supply is connected correctly.
4. If the water pressure is too low, please replace the earplugs with new ones to prevent blockage or damage, or select a stronger mode as needed.

Product Introduction

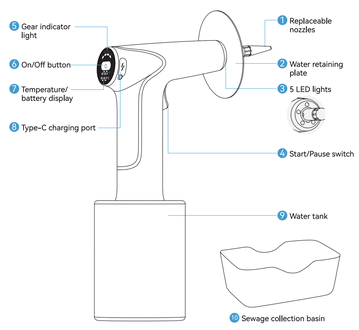Function Introduction

① Replaceable nozzle

The product is equipped with 1, 3, and 5-hole nozzles. You can choose to replace different nozzles according to actual conditions when using it. The more nozzle holes, the smaller the spray force

② Water baffle

When in use, you can install a water baffle to prevent water from spraying onto your body during use

③ 5 LED lights

When the product is working, the LED lights will automatically light up to illuminate the areas where work is needed

④ Start/pause switch

Press once to start, press again to pause, and use it in a cycle

⑤ Gear indicator light

The gears are divided into L, soft impact, M gear, medium impact, P gear, pulse mode, and H gear, strong impact. The corresponding gear light flashes when the work is suspended, and it is always on when working

⑥ On/off button

Press long to turn on/off, short press to switch the impulse gear, double-click to switch the power display

Switch between temperature units "°F and °C"

Long press the water switch (Figure 1), then short press the power switch (Figure 2) to switch the temperature unit.

⑦ Temperature and power display screen

The water temperature display corresponds to the "**°F °C**" icon on the side, and the power display corresponds to the "100" icon on the side. The "100" icon flashes when charging. The fully charged icon stops flashing and the battery level is displayed as "100".

⑧ Type-C charging port

Use a 5V-1A charger to charge the product, and use the Type-C interface for the charging cable

⑨ Clean water tank

Please fill it with warm water that adapts to body temperature when using it, and avoid using water that is too low or too clear, which may cause harm to the human body. The water tank capacity is 500 ml

⑩ Sewage collection tank

When using it, you can put the sewage collection tank next to your ear to collect sewage

Before use

1.It is normal when the product(water bottle, water tube, etc.)has some residual distilled water that used for testing performance.
2.It is recommended to use 1st gear mode for first time use,choose the other modes based on your needs.
3.It is recommended to soften heavy or hard earwax before irrigation.

Cleaning and maintenance

1.Take out the water tank, rinse the inside and outside of the water tank with water, and then let it air dry.

2. Replacement: After the nozzle has been used for 90 days, it needs to be replaced with a new one.
3. Storage: Wipe the device with a clean damp cloth. Make sure the device is dry before storing in a cool and dry place.

Notice:

(1) Avoid submerging the device in water to prevent damage.
(2) Avoid using harsh chemicals: Use dishwashing liquid or warm water for cleaning rather than harsh chemicals or solvents.

Environmental Protection

1. Proper disposal: When disposing of the flusher or its components, follow local e-waste regulations and guidelines and avoid disposing of it with normal household waste.
2. Battery recycling: Recycle rechargeable batteries properly. Check local facilities or collection points for battery recycling.
3. Energy efficiency: Save energy by turning off the flusher when not in use and unplugging the charger when the battery is fully charged.

Charging Instructions

*** Please use the Type-c charging cable that comes with the product to charge your device.**

1. Turn off the device before charging and make sure it is not in use.
2. Keep the device and charger away from water to prevent any electrical hazards.
3. The LCD display has a built-in power display for easy viewing of the power level.
4. After fully charged, disconnect the device from the charger.

Customer Support

If you have any questions about the performance of your product, or if you encounter a situation where accessories are missing or damaged, please contact our Customer Support.

Product Parameters

Product name: Powered Ear Wax Remover	Input: Type-C DC 5V/1A
Charging time: 3 hours	Water tank capacity: 500ml

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