

尺寸:85*125mm

BE011

User Manual for Automatic Ear Cleaner

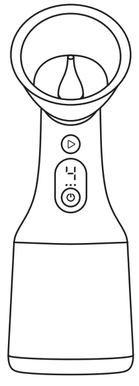

1. Appearance Overview

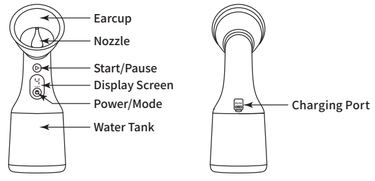

2. Accessories List




Main Unit*1 | Nozzles*4 | Charging Cable*1 | User Manual*1

3. Specifications

Product No: BE011	Product Size: 265*90*90mm
N.W.: 330g	Product Noise: ＜65dB
Rated Voltage: 3.7V	Water Tank Capacity: 200ml
Input Voltage: 5V⎓1A	Battery Capacity: 1000mAh

4. Function Description

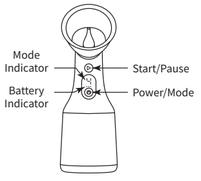


- **1st Level: Gentle Mode**
  (It is recommended to use the gentle mode for the first use)
- **2nd Level: Strong Mode**
- **3rd Level: Pulse Mode**
  (Pulse mode water spraying is intermittent.)
- **4th Level: Custom Mode**

In Custom Mode: Press and hold the "Start/Pause" button, and the level number will start flashing. At this point, the washing intensity will vary from low to high. Release the button when you feel the intensity is suitable. At this point, the 4th level will automatically remember this intensity, and later when you switch to the 4th level, it will be set to this custom washing intensity to accommodate different needs.

Battery: Full charge is indicated by three bars. When the three lights flash, it means the battery is low, and it should be charged promptly. When fully charged, all three bars will be lit.

5. Operating Instructions

*For first-time use, it is recommended to choose the gentle mode.
*It is recommended to soften earwax with a doctor-approved ear drop solution before washing.
*The device will automatically power off if not used within 10 minutes after startup.

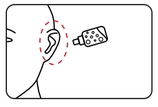


1. Drip cleaning solution into ear canal to soften earwax(purchase under doctor' s guidance)

Note: If the earwax is dry & hard, better gently loosen it with ear scoop first.




2. Install the nozzle, recommended to clean and disinfect the nozzle with 75% alcohol beforehand.

Note: When multiple people use the same machine, avoid cross-infection.

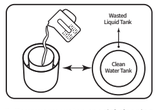

3. Pour warm water, special cleaning salt water or other solutions recommended by doctor into the clean water tank.

Note: Recommended to use warm water at 37°C±2°C or 98°F ±4°F

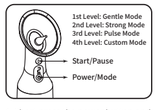

4. Short press the power button and select the appropriate mode.



5. Aim the nozzle at the ear canal and press the "Start/Pause" button to begin flushing. You can press the button at any time to stop during the flushing process.



6. After flushing, use a cotton swab or drying device to remove excess water from the canal to keep it dry.

6. Power On/Off, Pause, and Charging

Power On	Short press the Power/Mode Button to turn on the device. The indicator light will remain on, and the device will enter standby mode.
Pause	When the product is in operation, short press the Start/Pause button to pause. Press the Start/Pause button again to resume operation.
Power Off	Long press the Power/Mode Button for 2 seconds to turn off the device. The indicator light will go out.
Low Battery	When the three lights flash, it indicates low battery.
Charging	During charging, the charging indicator light: white light gradually flashing.
Charging Complete	Once charging is complete, the charging indicator light: white light stays on.

7. Precautions

Please read all instructions, warnings, and tips carefully before use. Improper use may cause personal injury and property damage.

1. Do not use this product if you have an ear infection, tympanic membrane perforation, or any ear injury.
2. If you experience pain during use, immediately stop using the product.
3. This product is strictly prohibited for use by children under 8 years old. Minors over 8 should use it under adult supervision.
4. This product is not a toy. Keep it out of the reach of children.
5. The product contains a built-in battery. Keep it away from flames, moisture, and high-temperature environments, as it may explode if overheated.
6. When not in use, ensure the power plug is unplugged.
7. Do not use water or cleaning solutions below body temperature, as they may cause dizziness or other discomforts.
8. Do not disassemble this product on your own.
9. Dispose of the product and its parts according to local electronic product disposal methods once their lifespan ends to protect the environment.
10. Do not use the product while it is charging.

8. Maintenance and Care

1. Cleaning: After each use, empty the remaining liquid from the tank and wipe clean with a towel. Store in a cool, ventilated place.
2. The main device can be wiped with a wet cloth; do not rinse directly with water.
3. The nozzle can be disinfected with 75% alcohol. Do not use high-temperature or high-pressure disinfection.
4. Avoid cross-infection when multiple people use the same device.

9. Troubleshooting

Problem	Cause	Solution
No water output or insufficient pressure	Low battery	Charge promptly
	Blocked filter head	Clean or replace the filter head
Cracked nozzle	Severely worn nozzle	Replace with a new nozzle
Device won't start	Motor not running	Ensure the product is sufficiently charged or contact our after-sales department and your dealer.
Other issues	Contact our after-sales department or your local dealer.	

CE FC RoHS Made In China