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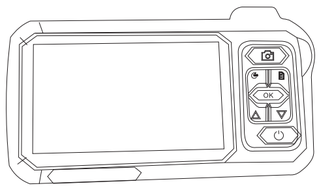
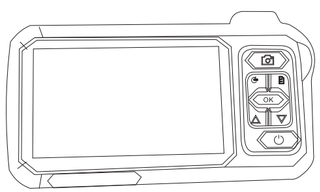
Model：T29


Warning

●You must read the instruction manual carefully before using.
●You must put the protective cover on to the otoscope probe before using, then sanitize all the areas that contact the skin.
●Keep out of reach of children. Children are not allowed to use this product alone, they must be used with the help of adults.
●Please take great care when using ear to clean earwax. DO NOT insert them too far in to the ear canal.
●Prevent someone from touching your hand when use the device in case of injury
●Stop using this product immediately if you experience any discomfort
●This product can not replace a doctor or medical providers diagnosis or treatment. Contact your doctor immediately if you experience any discomfort.

Battery Description

●Please fully charged for the first time use.
●Please use a 5V household charger that complies with safety regulations to charge this device.
●This device can work while charging
●If you do not use the device for a long time,please ensure that it is charged at least once every 3 months to prevent irreparable damage caused by excessive battery discharge.

Features

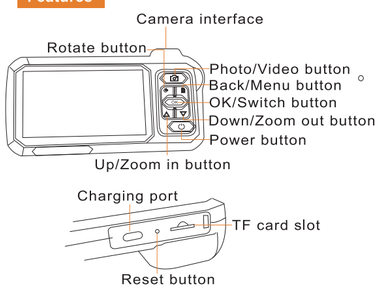Instructions for use

Connect the camera probe to the monitor by correct port

Extended press to the power button to turn on or off the device until see the LOGO image.

Tap the memory card to insert or eject.TF card status :T F card inserted/no card

The stand can be stored/-expanded, easy to use, free hands

Attention

●Power off device before inse rt TF card.
●Please place the TF card correctly.
●Do not use excessive force when inserting or flipping to avoid damaging the slot structure.

Language Setting

The default language is English, the system supports English, Chinese, German, French, Spanish, Italian, Russian and Japanese.
Setting path:
● Press the “” button to enter the setting page
●Select second item (Language)on setting page.
●Press up and down keys to select the language you need
●Press OK to save your choice.

Specifications

About the screen
Screen type: IPS display
Screen resolution(Pixel)： 854*480
photo resolution(Pixel) 1920*1080
Video Recording resolution(Pixe l) 1920*1080
Port: Type-C charging port, TF memory card slot

Camera
Diameter： 3.9mm
Camera resolution: 2 million pixels
Viewing angle:70°
Depth of Field: I5-20mm


Others
Operating Temperature:32°Fto I I3°F(0°C to 45°C)
Power Source:3000mAh Lithium Battery
Max battery life:4 ~5hours
Charging time:2~3 hours

FAQ

TROUBLSHOOTING GUIDE
Q: Why is the image so blurry? It used to be super clear.
AI: The No.I reason that cause those blurry images are the ear wax being stuck to the tip of the camera.Please use a cloth,alcohol,pad,or a Q-tip to clean it.You can apply a little pressure to make sure you get all the ear wax/dirt off the tip of the camera.
A2: The ear canal is not perfectly straight, gently adjust the position of the otoscope to view your eardrum clearly

Q: What should I do when images appear frozen or says "No camera"?
A: If the image is frozen,please REBOOT and RESTART the monitor. Finally, try to find "reset"and start it.Important, ensure the camera interface is fully inserted to the monitor.
Q: What is focal range?
A2: Please keep a focal length of about 1cm to achieve the best quality image(placing the scope too far from your focal point may result in blurry images)

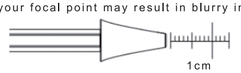Q: I bought it a few months ago and it worked well,but now I can't turn it on, looks like it does not work anymore.
A: Please make sure the battery is fully charged before you use it.If you cannot turn on the scope,or it turns off after a few minutes of working, it can be the sign of low battery.

⚠ Warning! Read all safety warnings and instructions before using this product. Failure to follow the warnings and instructionns may result in product damage.

MADE IN CHINA